



ABBOTT LABORATORIES HOSPITAL PRODUCTS DIVISION-MAN	
Figure 5-21. Power Supply (Continued) Schematic	
DRAWING NO.	REV. T
394-0206-010	SHEET 2 OF 2

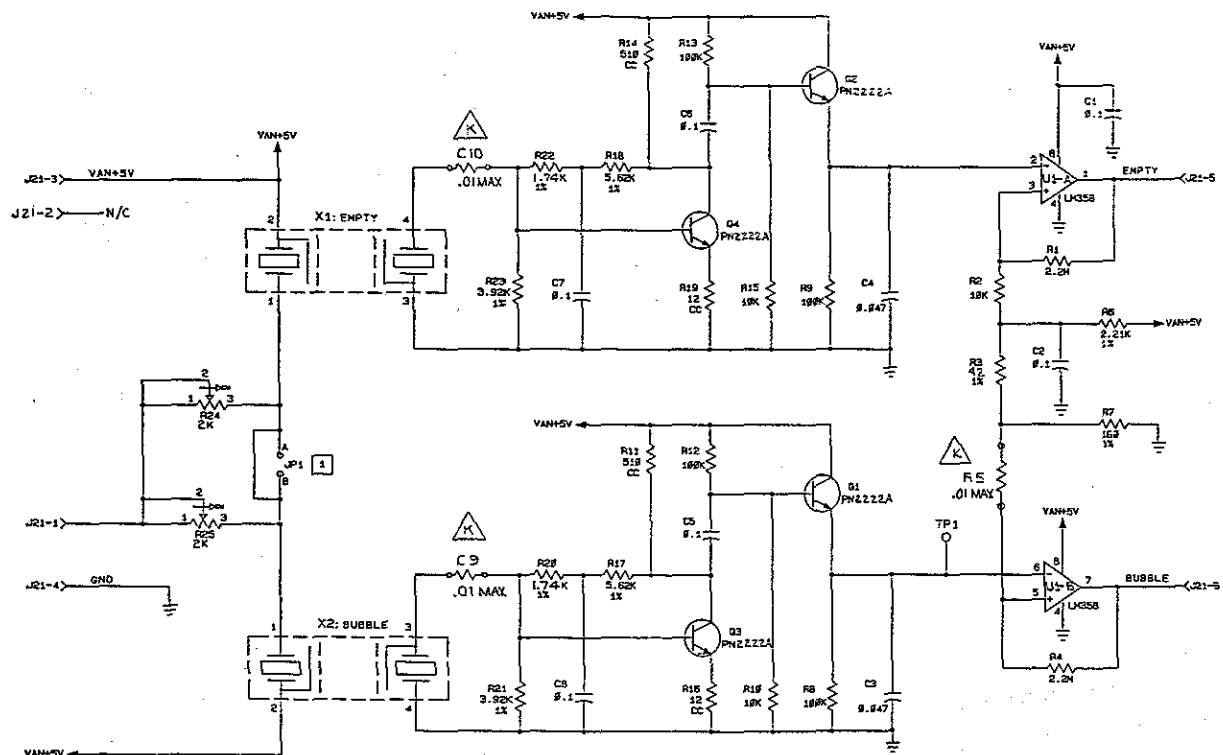

NOTES:

1. CUT TRACE BETWEEN JP1-A AND JP1-B IF R25 IS USED.
2. RESISTORS IN OHMS 5%, 1/4W UNLESS NOTED.
3. CAPACITORS IN MICROFARADS 25% UNLESS NOTED.

ABBOTT LABORATORIES HOSPITAL PRODUCTS DIVISION-MP	
Figure 9-22. Bubble Sensor (International) Schematic	
DRAWING NO.	REV. K
219-02222-010	SHEET 1 OF 1

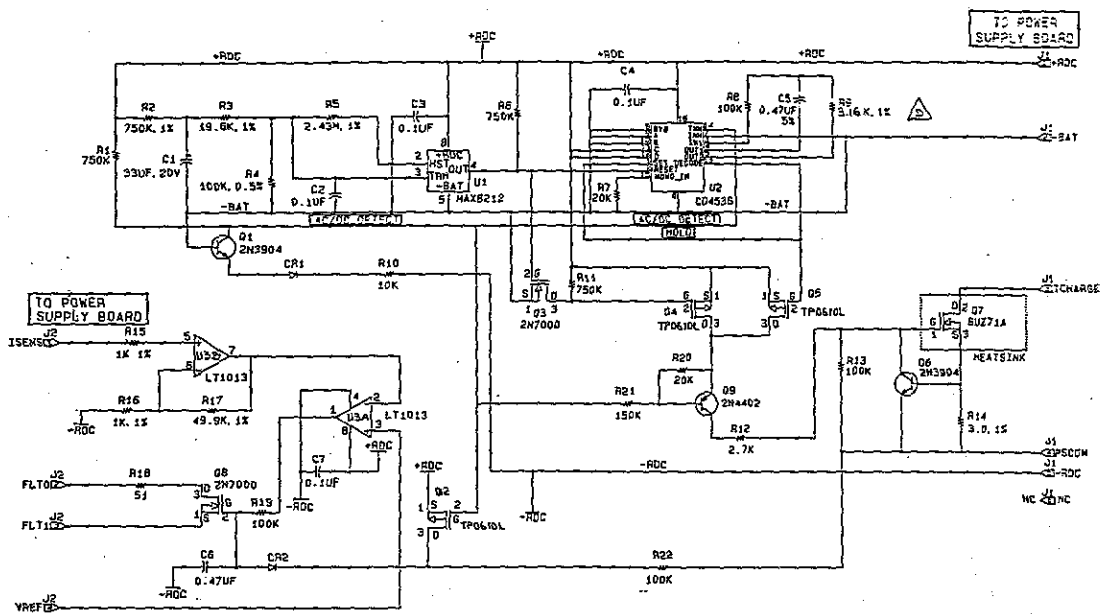

ASBOTT LABORATORIES HOSPITAL PRODUCTS DIVISION-MY	
Figure 1-23. Constant Bias Charger Schematic	
DRAWING NO. 94-2844-000	REV. D SHEET 1 OF 1

## Section 10

# INDEX

---

+5 VDC Supply Control, 4-31

-5 Volt Generation, 4-23

## A

A/D Conversion

1.5 Series, 4-23

1.6 Series, 4-23

AC (Mains) Line and Battery Power

Indication, 4-30

AC (Mains) Operation, 4-36

AC (Mains) Power Cord

Plug Replacement, 7-8

Replacement, IEC Qualified, 7-7

Replacement, UL Qualified, 7-6

AC Power Interface, 4-17

AC Receptacle Assembly Replacement,  
7-39

Accessories 1.6 Series, 7-4

Acronyms and Abbreviations, 1-4 - 1-5

Adjustment

Distal Pressure Sensor Electrical, 7-44

LCD Screen Contrast, 7-17

Air-in-Line Detection, 4-13

Distal, 4-13

in Cassette, 4-39

Proximal, 4-13

Air-in-Line Test

1.5 Series, 5-8

1.6 Series, 5-17

Alarm

and Malfunction Codes, 6-20

Conditions, 4-6

History Error Codes, 4-12

History, Obtaining 1.5 Series, 6-19

History, Obtaining 1.6 Series, 6-19

Messages, 6-2

Alarm Codes, 6-20

and Corrective Actions, 6-2

Display, Audible, and Touchswitch, 6-20

Infusion Pumping Mechanism, 6-20

Microprocessor or System, 6-20

Miscellaneous, 6-20

Alarm, Audible

Backup, 4-32

Control, 4-32

Self Test, 4-33

Alarms, Audible, 6-1

AM and Threshold Detectors, 4-42

Analog Multiplexer, 4-44

Artifacts, 1-7

Assistance, Technical, 6-1

## B

Bag Hanger Replacement, 7-14

Battery

Charger Operation, 4-9

Depth of Discharge, 4-7

Operation, 4-8

Operational Requirements, 4-8

Overview, 4-7

Power Indication, AC (Mains) Line, 4-30

Recharge, 4-8

Voltage Detection, 4-34

Battery Charger PWA, 4-35

Functional Block Diagram, 4-36

Interface, 4-17

Replacement 1.6 Series, 7-25, 7-27

Schematic 1.6 Series, 9-21

Battery Charger Test

1.5 Series, 5-8

1.6 Series, 5-17

Battery Pack

Charging, 4-33

Replacement, 7-5

Bubble Sensor Location Test

1.5 Series, 5-6

1.6 Series, 5-15

Bubble Sensor PWA, 4-40

Functional Block Diagram, 4-41

Schematic, 9-43

Bubble Sensor Schematic

(International), 9-65

Bubble Sensor Tips Removed

Infusion System Cassettes with

1.5 Series, 5-8

Infusion System Cassettes with

1.6 Series, 5-17

## C

Cassette, 4-47

Door Shaft Replacement, 7-34

Door Subassembly, 4-49

Dry 1.5 Series, 5-7

## INDEX

- Fluid Path in the, 4-2
- Installation Recognition, 4-38
- Leak Tests 1.5 Series, 4-3
- Leak Tests 1.6 Series, 4-4
- Loading Leaf Spring Replacement, 7-29
- Major Elements of the Dual-Channel, 4-42
- Character LEDs and Drivers, 4-44
- Cleaning Solutions, 5-2
- Cleaning the Infusion System, 5-2
- Clock
  - MPU and, 4-21
  - Second Source, 4-22
- Clutch
  - Housing Replacement, 7-14
  - Spring Replacement, 7-15
- Concurrent Delivery Test
  - 1.5 Series, 5-9
- Conventions, 1-4
- Core Failure State, 4-10
- Cotter Ring Replacement, 7-13
- Critical Data Retention, 4-11
- Current Boost Charger Schematic, 9-67
- Current Sensing Amplifier, 4-44
- Custom I/O IC, 4-24
- Custom Logic IC, 4-23

## D

- Data Retention, 4-11
  - Critical, 4-11
- DataPort
  - Accessory Cable Replacement, 7-40
  - Accessory Cable Schematics, 7-41
  - Interface, 4-18
  - Option, 4-26
  - Troubleshooting, 6-21
- DataPort Communication Test
  - 1.6 Series, 5-20
- DC Operation, 4-37
- DC Power Interface, 4-17
- DC Power Supply
  - Unregulated, 4-30
- Delivery Accuracy Test
  - 1.5 Series, 5-9
  - 1.6 Series, 5-18
- Delivery Mode
  - Selection, 4-26
  - Setting, 1-10 - 1-11
- Detectors, Monitors and, 4-12
- DIP Switch
  - Cover Replacement, 7-11
  - Settings for Each Delivery Mode, 1-11
- Display, Audible, and Touchswitch Alarm Codes, 6-20
- Distal Occlusion, 4-12

- Distal Pressure Sensor Electrical Adjustment, 7-44
- Door
  - Assembly Parts Replacement, 7-33
  - Assembly Replacement, 7-30 - 7-31
  - Base Casting Replacement, 7-35
  - Cassette, Shaft Replacement, 7-34
  - Cover Replacement, 7-32
  - Handle Replacement, 7-27
  - Handle, Mechanism and Leaf Spring Replacement, 7-28
  - Hinge Loading Spring Replacement, 7-33
  - Mechanism Leaf Spring Replacement, 7-28
  - Shield Replacement, 7-45
- Drawings, 9-1
- Dry Cassette, 5-7

## E

- Electrical Adjustment
  - Distal Pressure Sensor, 7-44
- Electrical Safety Test
  - 1.5 Series, 5-12
  - 1.6 Series, 5-20
- EMI Shield Replacement, 7-16
- Empty Container Test
  - 1.5 Series, 5-7
  - 1.6 Series, 5-15
- End of Performance Verification Test
  - 1.5 Series, 5-12
  - 1.6 Series, 5-21
- EPROM, RAM and Memory Protection
  - 1.5 Series, 4-22
  - 1.6 Series, 4-22
- Equipment Required
  - 1.6 Series, 5-13
  - 1.5 Series, 5-4
- Error Codes
  - Alarm History, 4-12
- Exiting from Failure State, 4-11

## F

- Failure State
  - Core, 4-10
  - Exiting from, 4-11
  - Peripheral, 4-10
- Features, Series Specific, 1-10
- Flow Detector 1.6 Series, 4-14 - 4-15, 4-26
  - Connected During Reset, 4-14
  - Disconnected While System is Operating, 4-15
  - Interface, 4-19
  - PWA, 4-44

- PWAs Schematic Diagram, 4-45
- Replacement, 7-12
- Fluid Path in the Cassette, 4-2
- Form
  - PVT Data 1.5 Series, 5-23
  - PVT Data 1.6 Series, 5-25
- Front and Rear Covers
  - Replacement, 7-16
  - Separating, 7-15
- Front Panel Assembly Replacement, 7-17
- Front Panel Touchswitches, 4-18
- Fuse and Fuse Drawer Replacement, 7-8

## G

- Gauge Dial Indicator, 5-6, 5-15

## H

- Handtools, Standard, 7-3

## I

- I/O Circuitry
  - Miscellaneous 1.5 Series, 4-26
  - Miscellaneous 1.6 Series, 4-27
- I/O Port Panel Replacement, 7-37
- I/O PWA, 4-23
  - Replacement, 7-21
  - Schematic 1.5 Series, 9-25
  - Schematic 1.5 Series, International), 9-57
  - Schematic 1.6 Series with DataPort), 9-29
  - Schematic 1.6 Series, 9-27
  - Schematic 1.6 Series, International), 9-59
  - with DataPort Option Functional Block Diagram, 4-24
- Infusion Mode Configurations, 1-2
- Infusion System
  - Cleaning the, 5-2
  - Inspecting the, 5-1
  - Sanitizing the, 5-3
- Infusion System Cassettes
  - with Bubble Sensor Tips Removed
    - 1.5 Series, 5-8
  - with Bubble Sensor Tips Removed
    - 1.6 Series, 5-17
- Inspecting the Infusion System, 5-1
- Inspection, 1-8
  - 1.5 Series, 5-5
  - 1.6 Series, 5-13
- Instrument Installation Procedure, 1-7
  - Inspection, 1-8
  - Self Test, 1-8

- Unpacking, 1-8
- Interconnect
  - Schematic 1.5 Series, 9-17
  - Schematic 1.6 Series, 9-19
- Interface
  - AC Power, 4-17
  - Battery Charger PWA, 4-17
  - Block Diagram, 4-16
  - DataPort, 4-18
  - DC Power, 4-17
  - Description, 4-15
  - Display, 4-19
  - Flow Detector 1.6 Series, 4-19
  - Main and I/O, 4-19
  - Motors and Valve, 4-18
  - Nurse-Call, 4-18
  - Power, 4-17
  - Power Supply PWA, 4-17
  - Pressure Sensor, 4-19
  - Sensor, 4-19
  - User, 4-18
- Introduction to Product, 1-1
- IPB for the Infusion System, 9-2, 9-7

## J

- Junction Box Assembly (DataPort Option), 4-46
- Junction Box PWA
  - Functional Block Diagram, 4-46
  - Schematic 1.6 Series, 9-51
- Junction Box Replacement
  - 1.6 Series with DataPort), 7-40

## L

- LCD Assembly, 4-45
  - Replacement, 7-20
- LCD Screen Contrast Adjustment, 7-17
- Leak Tests, Cassette
  - 1.5 Series, 4-3
  - 1.6 Series, 4-4
- LED Display PWA, 4-43
  - Functional Block Diagram, 4-43
  - Replacement, 7-19
  - Schematic, 9-23

## M

- Main and I/O Interface, 4-19
- Main Program Software Functions, 4-5
- Main PWA, 4-20
  - and I/O PWA Replacement, 7-21
- Functional Block Diagram, 4-21
- Replacement, 7-21
- Schematic 1.5 Series, 9-31

## INDEX

- Schematic 1.5 Series, International), 9-53
- Schematic 1.6 Series, 9-33
- Schematic 1.6 Series, International), 9-55
- Maintenance
  - Preventive, 5-1
  - Tests, 5-1
- Major Elements of the Dual-Channel Cassette, 4-42
- Malfunction
  - and Alarm Codes, 6-20
  - Detection, 4-10, 4-14
- Materials, 7-4
- Materials Required
  - 1.5 Series, 5-4
  - 1.6 Series, 5-13
- Mechanical Elements Behind Cassette Door, 5-3
- Mechanics Functional Description, 4-47
- Mechanism Assembly, 4-48
  - Cleaning, 7-42
  - Elements on the, 4-49
  - Lubrication, 7-42
  - Replacement 1.5 Series, 7-23 - 7-24
  - Replacement 1.6 Series, 7-25 - 7-26
- Mechanism Assembly Lubrication Points, 7-44
- Microprocessor or System Alarm Codes, 6-20
- Minipole Assembly Replacement, 7-13, 7-14
- Monitors and Detectors, 4-12
- Motor Drivers, 4-25
- Motor Power Supply Control, 4-31
- Motors and Valve Interface, 4-18
- Motors and Valves, 4-49
- MPU and Clock, 4-21

## N

- Nurse-Call
  - Cable Replacement, 7-13
  - Interface, 4-18
  - Relay Control 1.5 Series, 4-25
  - Relay Control 1.6 Series, 4-25
- Nurse-Call Test
  - 1.5 Series, 5-7
  - 1.6 Series, 5-15

## O

- Off Status, 4-1
- On Status, 4-2
- Operating Manual System, 3-1

- Operation, 4-6
  - Theory of, 4-1
- Operations
  - Sequence of, 4-1
- Optical Interrupters, 4-40
- Overvoltage Protection, 4-32

## P

- Parts
  - Replaceable, 9-2, 9-7
  - Replaceable, List, 7-1
  - Replaceable, Repairs and, 7-1
- Performance Verification Test
  - 1.5 Series, 5-4
  - 1.6 Series, 5-12
- Peripheral Failure State, 4-10
- Piezoelectric Alarm Assembly
  - Replacement, 7-38
- Plunger Drive Assembly, 4-50
- Plunger Motor Shaft and Nut
  - Lubrication, 7-44
- Pole Clamp
  - and EMI Shield Replacement, 7-10
  - Component Replacement, 7-9
  - Friction Plate Replacement, 7-11
  - Knob Replacement, 7-10
  - Shaft and Screw Replacement, 7-10
- Power Control, 4-30
- Power Interface, 4-17
- Power Supply
  - Schematic (International), 9-61
- Power Supply PWA, 4-28
  - Functional Block Diagram, 4-29
  - Interface, 4-17
  - Replacement 1.5 Series, 7-23
  - Replacement 1.6 Series, 7-25
  - Schematic 1.5 Series, 9-35
  - Schematic 1.6 Series, 9-39
- Preamplifier, 4-42
- Pressure Amplification, 4-39
- Pressure Sensing System, 4-12
- Pressure Sensor Interface, 4-19
- Pressure Sensor Test
  - 1.5 Series, 5-10
  - 1.6 Series, 5-19
  - Setup 1.5 Series, 5-11
- Preventive Maintenance, 5-1
- Primary/Secondary Valve Assembly, 4-50
- Proximal Occlusion
  - 1.5 Series, 4-13
  - 1.6 Series, 4-13
- PVT Data Form
  - 1.5 Series, 5-23
  - 1.6 Series, 5-25

## PWA

Battery Charger, 4-35  
Bubble Sensor, 4-40  
Flow Detector 1.6 Series, 4-44  
Functional Description, 4-20  
I/O, 4-23  
LED Display, 4-43  
Main, 4-20  
Power Supply, 4-28  
Sensor, 4-37

## Q

Qualification, User, 1-7

## R

Recirculating Set Test Setup 1.6 Series, 5-16  
Replaceable  
Parts and Repairs, 7-1  
Parts List, 7-1  
Replacement  
AC (Mains) Power Cord (IEC Qualified), 7-7  
AC (Mains) Power Cord (UL Qualified), 7-6  
AC (Mains) Power Cord Plug, 7-8  
AC Receptacle Assembly, 7-39  
Bag Hanger, 7-14  
Battery Charger PWA 1.6 Series, 7-25, 7-27  
Battery Pack, 7-5  
Cassette Door Shaft, 7-34  
Cassette Loading Leaf Spring, 7-29  
Clutch Housing, 7-14  
Clutch Spring, 7-15  
Cotter Ring, 7-13  
DataPort Accessory Cable 1.6 Series, 7-40  
DIP Switch Cover, 7-11  
Door Assembly, 7-30  
Door Assembly Parts, 7-33  
Door Base Casting, 7-35  
Door Cover, 7-32  
Door Handle, 7-27  
Door Handle and Door Mechanism Leaf Spring, 7-28  
Door Hinge Loading Spring, 7-33  
Door Mechanism Leaf Spring, 7-28  
Door Shield, 7-45  
EMI Shield, 7-16  
Flow Detector 1.6 Series, 7-12  
Front and Rear Cover, 7-16  
Front Panel Assembly, 7-17

Front Panel Assembly, Main PWA and I/O PWA, 7-19  
Fuse and Fuse Drawer, 7-8  
I/O Port Panel, 7-37  
I/O Port Panel and Piezoelectric Alarm, 7-38  
I/O PWA, 7-21  
Junction Box 1.6 Series With DataPort), 7-40  
LCD Assembly, 7-20  
LED Display PWA, 7-19  
Main PWA, 7-21  
Main PWA and I/O PWA, 7-21  
Mechanism Assembly  
1.5 Series, 7-23  
1.6 Series, 7-25  
Minipole Assembly, 7-13 - 7-14  
Nurse-Call Cable, 7-13  
Piezoelectric Alarm Assembly, 7-38  
Pole Clamp and EMI Shield, 7-10  
Pole Clamp Component, 7-9  
Pole Clamp Friction Plate, 7-11  
Pole Clamp Knob, 7-10  
Pole Clamp Shaft and Pole Clamp Screw, 7-10  
Power Supply PWA  
1.5 Series, 7-23  
1.6 Series, 7-25  
Procedures, 7-3  
Rubber Foot Pad and Foot Insert, 7-12  
Touchswitch Panel, 7-20  
Velcro Strap and Retaining Plate, 7-9  
Required  
Materials, 7-3  
Tools, 7-3  
Required 1.5 Series  
Equipment, 5-4  
Materials, 5-4  
Required 1.6 Series  
Equipment, 5-13  
Materials, 5-13  
Rubber Foot Pad and Foot Insert  
Replacement, 7-12  
  

## S

Sanitizing the Infusion System, 5-3  
Schematic  
Battery Charger PWA 1.6 Series, 9-21  
Bubble Sensor (International), 9-65  
Bubble Sensor PWA, 9-43  
Current Boost Charger, 9-67  
DataPort Accessory Cable, 7-41  
I/O PWA 1.5 Series, 9-25  
1.5 Series, International, 9-57  
1.6 Series with DataPort, 9-29



- 1.6 Series, 9-27
- 1.6 Series, International), 9-59
- Interconnect 1.5 Series, 9-17
- Interconnect 1.6 Series, 9-19
- Junction Box PWA 1.6 Series, 9-51
- LED Display PWA, 9-23
- Main PWA 1.5 Series, 9-31
  - 1.5 Series, International, 9-53
  - 1.6 Series, 9-33
  - 1.6 Series, International, 9-55
- Power Supply (International), 9-61
- Power Supply PWA 1.5 Series, 9-35
- Power Supply PWA 1.6 Series, 9-39
- Sensor PWA 1.5 Series, 9-45
  - 1.6 Series, 9-47
- Scope of Manual, 1-3
- Second Clock Source, 4-22
- Self Test, 1-8, 4-3
- Sensor Interface, 4-19
- Sensor PWA, 4-37
  - Functional Block Diagram, 4-38
  - Schematic 1.5 Series, 9-45
  - 1.6 Series, 9-47
- Separating the Front and Rear Covers, 7-15
- Series Specific Features, 1-10
- Service Tests, 5-1
- Setting the Delivery Mode, 1-10 - 1-11
- Setup, 4-5
- Software
  - Main Program Functions, 4-5
  - Obtaining Version Number, 1-9
- Solutions
  - Cleaning, 5-2
- Specifications
  - Domestic Infusion System 1.5 Series, 8-1
  - Domestic Infusion System 1.6 Series, 8-4
  - International Infusion System
    - 1.5 Series, 8-7
    - 1.6 Series, 8-10
- Start-Up Test
  - 1.5 Series, 5-5
  - 1.6 Series, 5-14
- Supply Control
  - +5 VDC, 4-31
  - Motor Power, 4-31
  - VMEM, 4-31
- System Malfunction Detection, 4-10
- T
  - Technical Assistance, 6-1
  - Test
    - Air-in-Line 1.5 Series, 5-8
    - 1.6 Series, 5-17
    - Battery Charger 1.5 Series, 5-8
    - 1.6 Series, 5-17
    - Bubble Sensor Location 1.5 Series, 5-6
    - 1.6 Series, 5-15
    - Concurrent Delivery 1.5 Series, 5-9
    - DataPort Communication 1.6 Series, 5-20
    - Delivery Accuracy 1.5 Series, 5-9
    - 1.6 Series, 5-18
    - Electrical Safety 1.5 Series, 5-12
    - 1.6 Series, 5-20
    - Empty Container 1.5 Series, 5-7
    - 1.6 Series, 5-15
    - End of Performance Verification
      - 1.5 Series, 5-12
      - 1.6 Series, 5-21
    - Nurse-Call 1.5 Series, 5-7
    - 1.6 Series, 5-15
    - Performance Verification 1.5 Series, 5-4
    - 1.6 Series, 5-12
    - Pressure Sensor 1.5 Series, 5-10
    - 1.6 Series, 5-19
    - Self, 4-3
    - Start-Up 1.5 Series, 5-5
    - 1.6 Series, 5-14
    - Test Setup
      - Pressure Sensor 1.5 Series, 5-11
      - Recirculating Set 1.6 Series, 5-16
    - Tests
      - Maintenance, 5-1
      - Service, 5-1
      - Threshold Detectors
      - AM and, 4-42
      - Touchswitch Panel Replacement, 7-20
      - Touchswitches, Front Panel, 4-18
      - Troubleshooting, 6-1, 6-21
      - DataPort Systems, 6-21
      - Troubleshooting with the PVT
        - 1.5 Series, 6-23
        - 1.6 Series, 6-25
  - UVW
    - Unpacking, 1-8
    - Unregulated DC Power Supply, 4-30
    - User
      - Interface, 4-18
      - Qualification, 1-7
      - Valves and Motors, 4-49
      - Velcro Strap and Retaining Plate
        - Replacement, 7-9
      - VMEM Supply Control, 4-31
      - Warranty, 2-1

For technical assistance within the United States, contact:

Abbott Laboratories  
Technical Support Operations  
960 Linda Vista Avenue  
Mountain View, California 94043

CALL TOLL FREE: 1-800-241-4002

For technical assistance outside the United States, contact the local sales office.

**CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician or other licensed practitioner.**

## WARNING

**POSSIBLE EXPLOSION HAZARD EXISTS IF USED IN THE PRESENCE OF  
FLAMMABLE ANESTHETICS.**

LifeCare, DataPort, PlumSet, and Plum are registered trademarks of Abbott Laboratories. Super Edisonite is a registered trademark of S. M. Edison Chemical Co. Velcro is a registered trademark of Velcro USA, Inc. X-acto is a registered trademark of Hunt Mfg. Braycote is a registered trademark of Castrol, Inc. Vesphene IIse and Manu-Klenz are registered trademarks of Calgon Vestal Laboratories. Formula C is a trademark of Diversey Corporation. Bio-Tek is a registered trademark of Bio-Tek Instruments, Inc. Dynatech Nevada is a registered trademark of Dynatech Nevada, Inc. Fluke is a registered trademark of John Fluke Manufacturing Company, Inc. Electro-Wash is a registered trademark of Chemtronics Inc.

Covered by one or more of the following U.S. patents: 4,336,800; 4,453,932; 4,552,336

IEC 601-1 Classification: Class 1, Type CB Drip-proof Medical Equipment

## IPX1



Attention, consult accompanying documents.



LISTED  
MEDICAL  
EQUIPMENT  
873K



LR 36894  
® RISK CLASS 2G  
C 22.2 No. 125

Listing does not cover IV sets



NRTL/C  
CSA 22.2  
UL 544

The use of NRTL/C adjacent to the CSA mark indicates that the product has been certified by CSA to U.S. and Canadian standards. CSA has been accredited by the U.S. Occupational Safety and Health Administration (OSHA) as a Nationally Recognized Test Laboratories (NRTL).

Printed in USA

Abbott Laboratories □ North Chicago □ Illinois 60064 USA